



Topi-click Troche-click, 0.90mL Cavity with Scored Halves and Quarters, Pink

Description: The Topi-Click Troche-Click is an innovative way to provide an accurate dose to patients in a troche lozenge. It is easy-to-fill, increases pharmacy efficiency and improves patient experience with medications dispensed directly from the disposable mold design.

Features:

- Each 0.90 mL cavity allows for customizable dosing with scored halves and quarters
- Easy-to-eject with eject button for quick removal
- Heat resistant tray
- Blocks UV
- Large label surface on cover
- Locking mechanism designed to prevent accidental opening
- Multiple color options make it easy for patients to distinguish their medications
- Materials compliant with FDA CFR 21.177 code, which is the requirement for <USP 661>

Technical Specifications:

Capacity	30 x 0.90 mL
Material	Polypropylene
Label Size	
Cover	3.0" W x 2.0" H
Tray	1.5" W x 0.5" H

USA | Canada | Australia | International


medisca.com

Filling Instructions



STEP 1:
Place the Troche-CLICK™
tray on a clean, level
surface.

Fill the Troche-CLICK™ with
mixture. Filling techniques
may include: pouring
directly from a glass
beaker, transferring mixture
using a syringe or transfer
pipette.

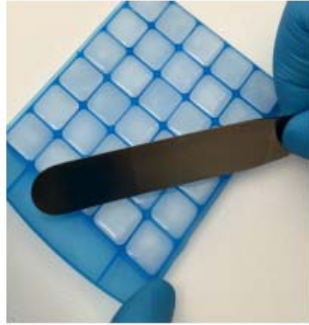

STEP 2:
Remove excess troche
material by scraping over the
top of the troche lozenges
with a clean spatula.

The PEG base shrinks
during cooling. Add
additional material to the
individual cavities, if required.



STEP 3:
After the troche material
has set, a heat gun can
be used to remelt the top
of the troches to provide
a polished appearance to
the troches.



STEP 4:
Troches lozenges can be
removed and then replaced
to conduct quality checks.

Related products:

- Topi-Click Troche-Click, 0.90mL Cavity with Scored Halves and Quarters, White (5494)
- Topi-Click Troche-Click, 0.90mL Cavity with Scored Halves and Quarters, Green (5493)
- Topi-Click Troche-Click, 0.90mL Cavity with Scored Halves and Quarters, Blue (5496)

USA | CANADA | AUSTRALIA PACIFIC | EUROPE
WWW.MEDISCA.COM

Please note that the information contained herein is for informational purposes only. The information is based on recognized compendia, manufacturers' instructions and/or other authoritative sources and texts. MEDISCA takes no responsibility for the validity or accuracy of this information. Customer MUST refer to USP, FDA, Health Canada, TGA and/or their State Board to ensure they comply with the specific regulations and ensure that appropriate procedures are in place.