



# Medi X-TABB™ Mold Cheat Sheet

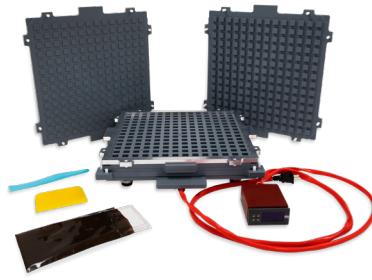

**Medi X-TABB™ Mold**  
Product No. 5650

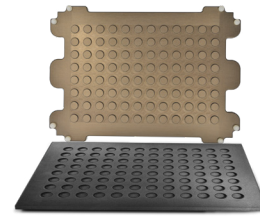

**ProFiller Rapid Dissolve Tablet Mold**  
Product No. 6707

<b>Appearance</b>	13 x 13 inches (L x W)	13 x 9.3 inches (L x W)
<b>Tablet Shape</b>	Square	Circular
<b>Cavity Volume</b>	1.024 mL	0.940 mL
<b>Material</b>	Non-stick, hard coated aluminum.	Teflon coated Hard anodized aluminum
<b>Cavities</b>	196	96
<b>Key Mold Components</b>	<ul style="list-style-type: none"> <li>• Cavity Plate</li> <li>• Tamper Plate</li> <li>• Quarter Score Plate</li> <li>• Base Heating Plate</li> </ul>	<ul style="list-style-type: none"> <li>• Cavity Plate</li> <li>• Tamper Plate</li> </ul>
<b>Baking</b>	Built-in base heating plate (does not require convection oven)	Requires Convection oven
<b>Baking Time</b>	110 °C (230 °F) for 15-20 mins (Medi X-TABB™ Base) 110 °C (230 °F) for 20-25 mins (Medi-RDT™ Base)	110 °C (230 °F) for 15-20 mins
<b>System Suitability</b>	This system increases compliance with USP 800 and makes the process efficient when handling HD APIs as entire compounding process, from filling to heating, can be performed under the hood.	Requires precautions with Hazardous APIs as while the mold can be fitted under the hood for filling and compression, the baking process requires the use of a convection oven that is typically outside of the hood.
<b>Blister Packing</b>	(5632) Medi X-TABB™ Blister Pack, 30 cavities, transparent PVC, requires sleeve or amber zip-lock bag for light resistant API's	(7219) Medi-RDT Blister Pack , 20 cavities, amber colored, light-resistant, perforated individual cavities with tamper evident seal

## Sample Formulations

<b>Medi X-TABB™ Mold + Medi X-TABB™ Base</b>	F 009 316	Progesterone 150 mg Oral Disintegrating Tablets
	F 009 317	Testosterone 20 mg Oral Disintegrating Tablets
	F 009 318	Estradiol 0.4 mg, Estriol 1.6 mg Oral Disintegrating Tablets
	F 009 319	Dehydroepiandrosterone 50 mg Oral Disintegrating Tablets
<b>Medi X-TABB™ Mold + Medi RDT™ Base</b>	F 009 330	Progesterone 200 mg Oral Rapid-Dissolve Tablets
	F 009 331	Melatonin 5 mg Oral Rapid-Dissolve Tablets
<b>Convection Oven + Medi RDT™ Base</b>	F 005 469v4	Progesterone 200 mg Oral Rapid-Dissolve Tablets
	F 005 472v4	Melatonin 5 mg Oral Rapid-Dissolve Tablets