


# Medi X-TABB™ Mold Sell Sheet

The Medi X-TABB™ Mold is a versatile and innovative system for preparing personalized tablets. Consisting of a calibrated mold with a built-in heating plate, this molding system allows pharmacies to compound oral dissolving tablets without the need of an oven! Thanks to its built-in heating element, the Medi X-TABB™ Mold can be used exclusively under a primary containment system making it ideal for hazardous drug compounding.

Built of durable and lightweight aluminum that is coated with the latest US FDA and EU compliant non-stick technology, the Medi X-TABB™ Mold provides an easy-to-use and easy-to-clean surface. The heating element will provide fast and even heating to ensure your tablets are uniformly baked while also providing the necessary temperature range to use several different bases. With the combination of the heating element and the aluminum material, the Medi X-TABB™ mold allows you to compound tablets in a quick and efficient manner to increase cycle times.

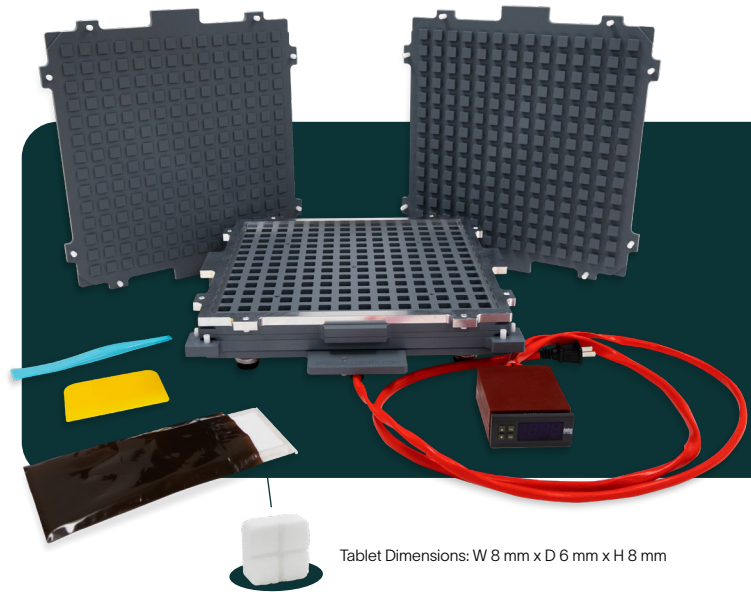

## Features

196 cavities with a volume of 1 mL

Tablets produced will be scored in quarters for added versatility

Fast and even heating (heats up to 60 °C in 30 seconds)

Celsius or Fahrenheit display options

Maximum heating temperature up to 260 °C (500 °F)

Ideal with Medi X-TABB™ Base

Also compatible with Medi-RDT™ Base

Approximate mold calibration for blank tablets:

- Medi X-TABB™ Base: 773 mg
- Medi-RDT™ Base: 724 mg

## Includes

Base Heating Plate with Digital Controller

Quarter Score Plate

Cavity Plate (196 cavities of 1 mL)

Tamper Plate

Scraper Spatula

Scraper/Pry Tool

Stainless Steel Thumb Screws

Instructional Manual

## Sample Formulations

<b>Medi X-TABB™ Mold + Medi X-TABB™ Base</b>	F 009 316	Progesterone 150 mg Oral Disintegrating Tablets
	F 009 317	Testosterone 20 mg Oral Disintegrating Tablets
	F 009 318	Estradiol 0.4 mg, Estriol 1.6 mg Oral Disintegrating Tablets
	F 009 319	Dehydroepiandrosterone 50 mg Oral Disintegrating Tablets
<b>Medi X-TABB™ Mold + Medi RDT™ Base</b>	F 009 330	Progesterone 200 mg Oral Rapid-Dissolve Tablets
	F 009 331	Melatonin 5 mg Oral Rapid-Dissolve Tablets