



Introducing

THE POWER OF PROMILL™

For The High Volume
Compounder




POWER • CALIBRATION • CLEANING

A MEDISCA Exclusive Innovation

THE NEW POWERFUL PROMILL™

MEDISCA's ProMill elevates ointment mills to a higher level of production. Its distinctive features include a capacity of up to 11 liters per hour, a self-calibration option and an innovative cleaning tool for its fixed rollers. This robust machine was designed specifically for compounders who need to mill large volumes of preparations in the most efficient manner.

ProMill is your new powerhouse mill.

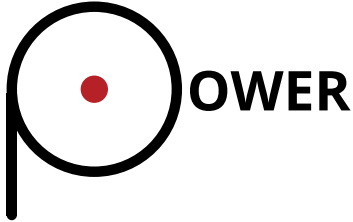


Features & Benefits

- High volume throughput, up to **11 liters per hour**
- Large **hopper of 1.75 liters** enables to process more material
- Digital control panel with variable speed control and counter time capabilities
- Solid and easy to clean **aluminum oxide rollers** for long term performance
- **Self-calibrating** gap distances
- Equipped with safety features: an emergency stop and engine overload protection

Excellent Return on Investment for High Volume Compounders

- Hourly Cost of Technician Labor and Benefits: \$25.00
- Assumption: Compounding 100 x 240 mL scripts = 24 L/day

Features	MEDISCA ProMill	Other Mill	ProMill Advantage
Roller Length	200 mm	< 150 mm	✓
Roller Diameter	50 mm	50 mm	=
Roller Material	Aluminum Oxide	Ceramic	✓
Max Output	11 L/h	7 L/h	✓
Hopper Capacity	1750 mL	700 mL	✓
Particle Fineness	< 15 µm	< 20 µm	✓
Interface	Digital	Analogue	✓

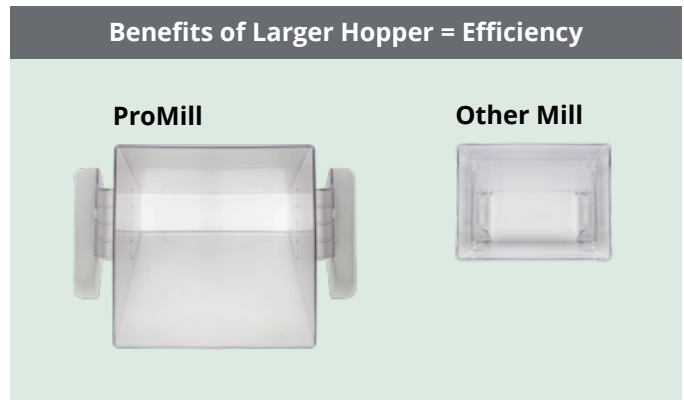


Model	Processing Rate	Processing Time (24 L daily)	Processing Time (260 work days per year)	\$/year (\$25 per Hour)
ProMill	11 L/h	2 hrs 15 min	585 hrs	\$14,625
Other Mill	7 L/h	3 hrs 30 min	910 hrs	\$22,750
LABOR SAVINGS			325 hrs	\$8,425
			RETURN ON INVESTMENT	9 months

Model	Processing Rate	Daily Processing Time	Number of 240 mL Scripts Processed
ProMill	11 L/h	4 hrs	183
Other Mill	7 L/h	4 hrs	116
INCREASED EFFICIENCY			67 Additional Scripts



TÜV SUD Certified for Your Safety

MEDISCA's ProMill has been certified by TÜV SUD for having successfully completed the testing against the standard IEC 61010-1: Safety Requirements for Electrical Equipment of Measurement, Control, and Laboratory Use. TÜV SUD is recognized by:

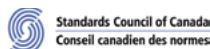

CALIBRATION

Take Power Into Your Own Hands!

In only **10 minutes**, you can calibrate the gaps should any misalignment occur. The calibration tool kit is supplied with each MEDISCA ProMill™ Starter Package, allowing you to control gap fineness and self-calibrate to avoid any downtime. **Increase efficiency by saving time and money.**

CLEANING

Ingenious Innovation!

Simplify and expedite your cleaning routine. MEDISCA's ProMill cleaning tool is an incredible, one of a kind cleaning device that facilitates the cleaning operation, and is **specifically designed for ProMill's fixed rollers.**

Cleaning Tool Microfibers provide additional convenience during the cleaning process.

It is as simple as Attach. Clean. Discard. Mill Again.



REAL TIME SAVER



Calibration Tool Kit

Product No.: 6380



Cleaning Tool Microfiber

Product No.: 6381

PK/25, PK/100, PK/500



Cleaning Tool

Product No.: 6382

Patent pending

STARTER PACKAGE INCLUDES

(Product No.: 6260-01)

Watch our New Video
at MEDISCA.com 

- 1. Hopper**
Product No.: 6385
- 2. Roller Guides**
Product No.: 6383
- 3. Calibration Tool Kit**
Product No.: 6380
- 4. Cleaning Tool**
Product No.: 6382
- 5. Microfiber Insert**
Product No.: 6381
- 6. Apron**
Product No.: 6384
- 7. Splash Tray**
Product No.: 6386



LIMITED TIME OFFER!



SAVE \$1000*

TRADE-IN any three roll ointment mill and upgrade to the **POWERFUL MEDISCA ProMill™**. Call your Account Executive today for more details.

* Cannot be combined with any other promotion or discount. Limited time offer. Conditions may apply. Applicable for US customers. Savings are based on catalog price.

 1.800.932.1039 WWW.MEDISCA.COM |  1.800.665.6334 WWW.MEDISCA.CA

FOLLOW US ON SOCIAL MEDIA    

 **MEDISCA**[®]
YOUR TRUSTED PARTNER IN COMPOUNDING