



# **QUICK JUST BECAME *FASTER!***



**WATCH OUR VIDEO**  
[www.medisca.ca/videos](http://www.medisca.ca/videos)

Now is the time to take your compounding practice to the next level by offering your patients MEDISCA's Rapid Dissolve Tablets. Our innovative RDT Kit accommodates your customers' needs.

## > WHAT ARE RDTs?


Rapid Dissolve Tablets are designed to dissolve quickly on the tongue without any water making it an ideal alternative dosage form for patients who have difficulty swallowing. In short, they offer quick delivery of medication, rapid dissolution, and greater absorption in the blood.

## > WHAT'S INCLUDED IN THE RDT KIT



1

- Medi-RDT Base
- The Medi-RDT Base is a Sucrose-Free, finely granulated powder that is compatible with a wide range of active ingredients. Our unique RDT formula provides significantly faster wetting and disintegration times compared to similar products on the market:
- Wetting time: 5 sec. (versus 24 sec.)
- Disintegrating time: 93 sec. (versus 287 sec.)
- Included in the kit: 100 g

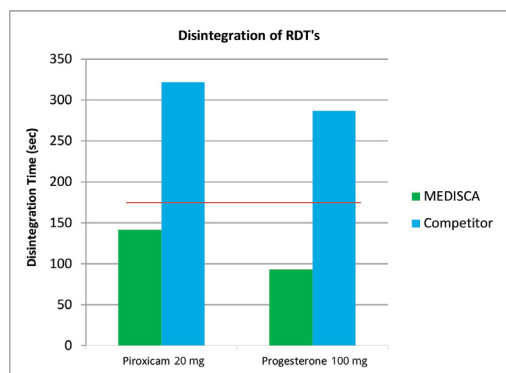


2

### PROFILLER RDT MOLD

The ProFiller RDT Mold is developed to provide quick, convenient and consistent medicated dosages. When used in conjunction with MEDISCA's Medi-RDT Base, it provides an easy and efficient way to dispense your RDTs with sufficient mechanical strength and quick disintegration in the mouth.

Available Mold:

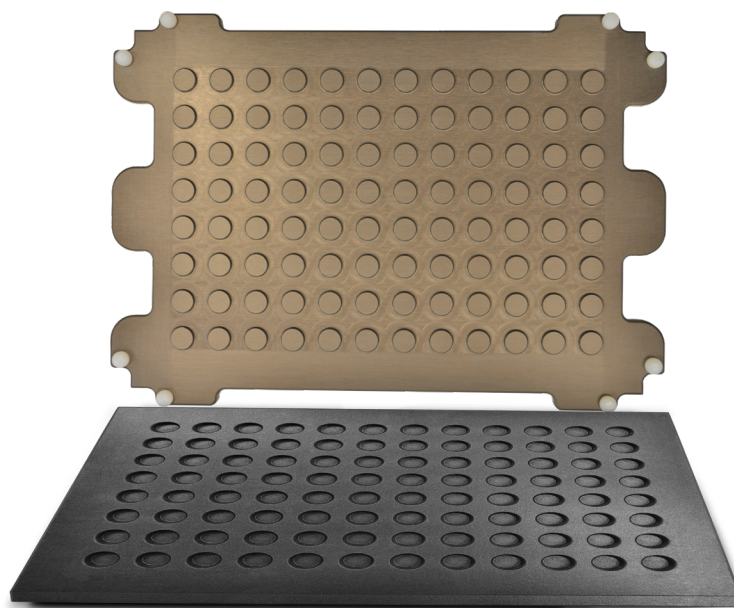

**6707-01: 1 PROFILLER RAPID DISSOLVE TABLET MOLD WITH 96 CAVITIES (0.94 ML/CAVITY)**

Included in the MEDI-RDT Kit (7218-01)

## > WHY CHOOSE MEDISCA's Medi-RDT BASE?

MEDISCA's Medi-RDT Sucrose-Free Base has been put to the test against competitors and has proven to have superior performance. This is one of the reasons for which more and more compounding pharmacists are choosing MEDISCA's Medi-RDT Base



3

### MEDI-RDT BLISTER PACKS WITH SLEEVES FOR 160 DOSES

The customized blister packs include UV-resistant amber-colored sheets for optimal protection of light-sensitive drugs. Cold seal tamper-evident foil laser labels for medication identification and sleeves are also included for easy dispensing.



4

### HOT HAND® PROTECTOR

With safety being a top priority, the Hot Hand protector mitt provides superior protection when handling hot or cold objects. The gripping surface is studded with suction cups for a secure no-slip grasp to safely remove the RDT Mold from the convection oven.



5

### CUSTOMIZED FORMULA

MEDISCA Network recognizes your need to get it right the first time which is why in order to ensure proper compounding method techniques for the Rapid Dissolve Tablets, MEDISCA Network will provide you with a customized sample formula of your choice to help get you started with confidence.

6

### 50 MESH SIEVE AND RECEIVER PAN

Another essential tool included is the 50 mesh sieve and pan which serve to control particle size distribution as well as the flowability of the blend. This is to ensure content uniformity of the solid dosage form.



# > GET STARTED TODAY!

Grow your practice and competitive edge with MEDISCA's RDT Kit. We want to set you up for success as quickly as possible with all the necessary tools to elevate your compounding practice and surpass your patients' expectations with our Rapid Dissolve Tablet Kit.

MEDISCA wants you and your patients to experience how **QUICK JUST BECAME FASTER!**



## MEDI-RDT KIT

<b>PRODUCT No.</b>	7218-01
<b>INCLUDES</b>	<ul style="list-style-type: none"> <li>• 1 x MEDI-RDT Base, 100 g</li> <li>• 1 x ProFiller Rapid Dissolve Tablet Mold, 96 Cavities x 0.94 mL</li> <li>• 1 x MEDI-RDT Blister Pack 20 Cavities w/Cold Seal &amp; Sleeves, 160 Doses</li> <li>• 1 x Hot Hand® Protector</li> <li>• 1 x Sieve &amp; Receiver Pan, 50 Mesh</li> <li>• 1 x Powder Scraper</li> </ul>

## REQUIRED ACCESSORIES

### 1 CONVECTION OVEN, FISHER SCIENTIFIC (Product No.: 6730-01)

- W 20.9" x D 22.2" x H 28.3", 482°F / 250°C, Digital Temperature, Forced Air

### 2 GC SERIES CONVECTION OVEN, QUINCY (Product No.: 8734-01)

- W 20.0" x D 14.3" x H 25.5", 450°F / 232°C, Analog Temperature, Gravity Convection



🇺🇸 1.800.932.1039 | 🇨🇦 1.800.665.6334 | 🇺🇰 1300 786 392 | 🌐 (01) 514.333.7811

WWW.MEDISCA.COM



DISCLAIMER: For use in prescription pharmacy compounding only, in accordance with applicable law. MEDISCA makes no warranty or representation regarding the use of this product or fitness for any particular purpose, and any such use is entirely the responsibility of the pharmacist in the exercise of his/her professional expertise and judgment. Accordingly, in no event will MEDISCA be liable for any damages from the use or misuse of the product. Purchaser's sole remedy in the event of product defect or other damage, shall be the refund of the purchase price of the product.